

Talkbox Project  
EECS 2311  
Alberto M., Rohan T., Jordan M.

## *TalkBox Configurator* User Manual


**T**alkbox is a highly customizable virtual lexigram interface that allows persons with speech disabilities to communicate. The configurator is rich with features to customize the appearance, functionality, feel of the audio buttons.

### 1 Installation Instructions

*Prerequisite:* **Java** must be installed. Click [here](#) to install.


To install the software, please visit the *Talkbox* project site at <https://github.com/Jakjm/Talkbox> and download and extract the zip folder. In the **Executables** folder, there are three core applications: The *Talkbox Configurator*, *Talkbox Simulator*, and *TBC Log*.

*Under the ‘Clone or download’ dropdown, click Download Zip.*

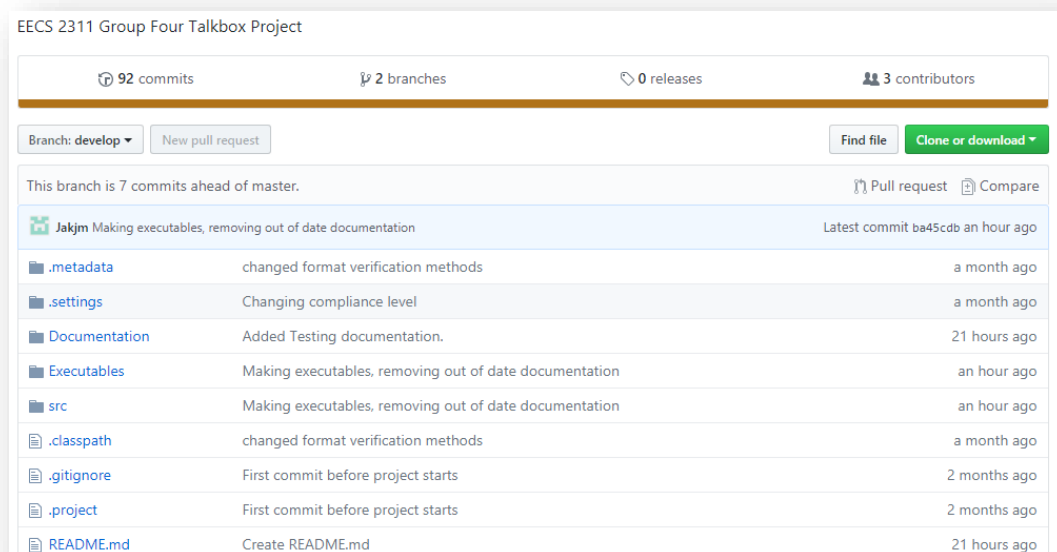

## 2 Using the Configurator

*In this configuration menu there are five buttons.*

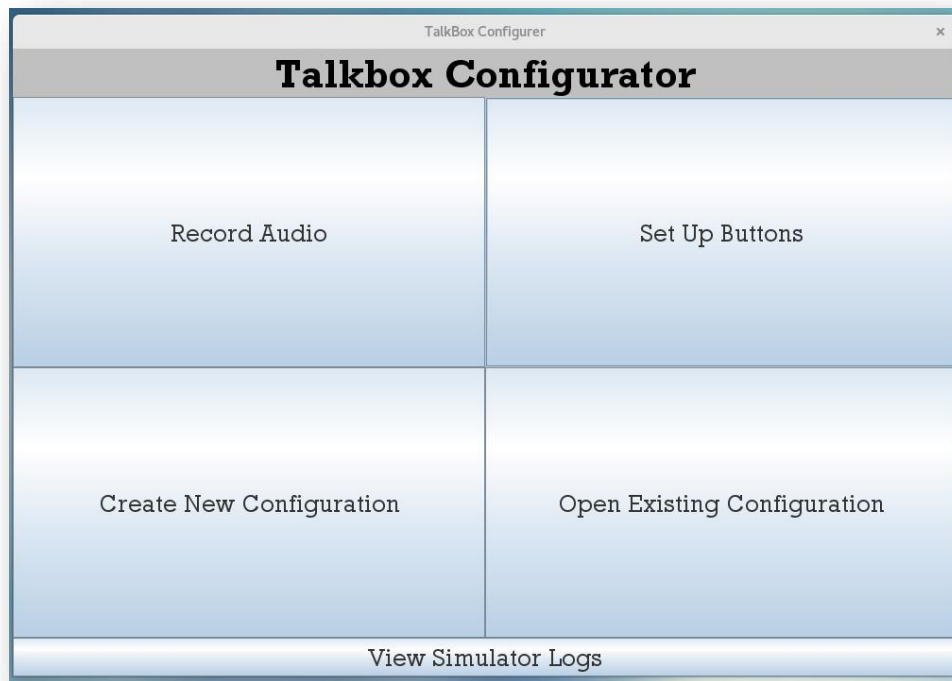

**Record Audio**—the purpose of this panel is to record phrases and words that can be used for the buttons. After recording, the user must save the audio file on their computer (but not the *Talkbox* directory— it will not recognize it). Alternatively, these sounds can be recorded within the *Set Up Buttons* panel.


The user can name the audio file to whatever they want. At this time only *WAVE* files can be used. To record audio, press *Record*. This is now the time to record anything that the user would like to record. When finished recording the user must press *Stop Recording*. After selecting *Stop Recording*, select a directory to save the recording in. When “OK” is selected, a window with a list of directories will pop up. The user should select where they would like to save the recording. This recording can then be played back when associated to a button. The below figure shows a screenshot of the audio recording panel.

Talkbox Project

EECS 2311

Alberto M., Rohan T., Jordan M.

*This is what the Recording audio panel looks like when selected*

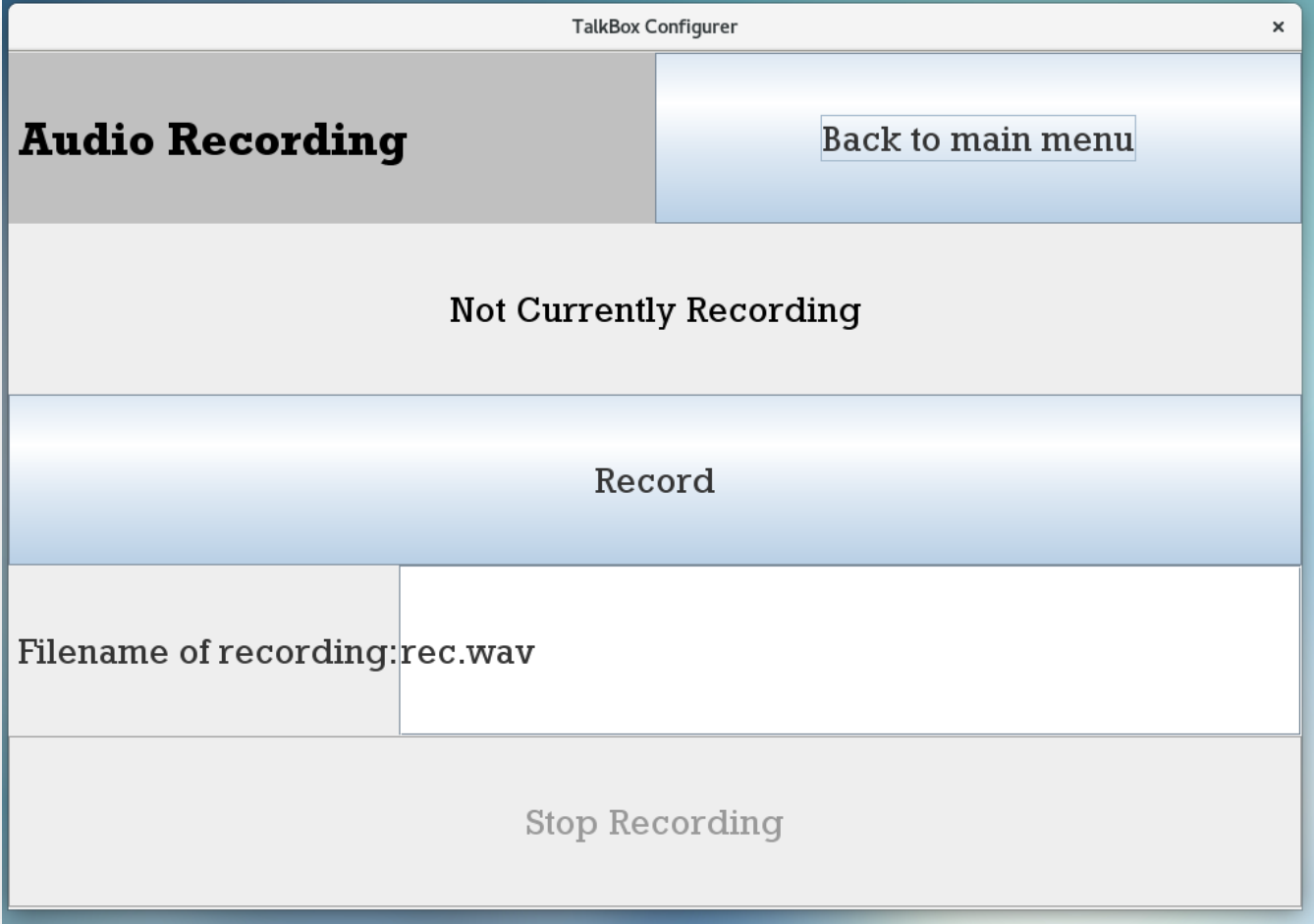


The screenshot shows a window titled "TalkBox Configurer" with a close button (x) in the top right corner. The window is divided into several sections. The top left section is a dark gray header with the text "Audio Recording" in bold white font. To its right is a light blue button labeled "Back to main menu". Below the header is a large light gray area with the text "Not Currently Recording" in bold black font. Underneath this is a wide light blue button labeled "Record". Below the "Record" button is a section with a light gray label "Filename of recording:" and a white text input field containing the text "rec.wav". At the bottom of the window is a wide light gray button labeled "Stop Recording".

TalkBox Configurer	
<b>Audio Recording</b>	Back to main menu
Not Currently Recording	
Record	
Filename of recording:	rec.wav
Stop Recording	

The **Create New Configuration** button allows the user to create a new *TalkboxData* directory. The **Open Existing Configuration** allows the user to select an existing *TalkboxData* configuration to allow further editing and customization.

The **Set Up Buttons** panel is not available until the user selects an existing *TalkboxData* directory or makes a new configuration directory. This is where one can setup the orientation of their buttons and can add or remove additional rows of buttons as one sees fit. In this panel the user can setup which audio clip goes to which button. They can use emojis, i.e. smile (😊) to allow the buttons to be more descriptive. The user can also name the button, change its color, or even add an image. When doing all these edits to the button, to ensure that this will be saved the user **MUST** press *Confirm Setup*.

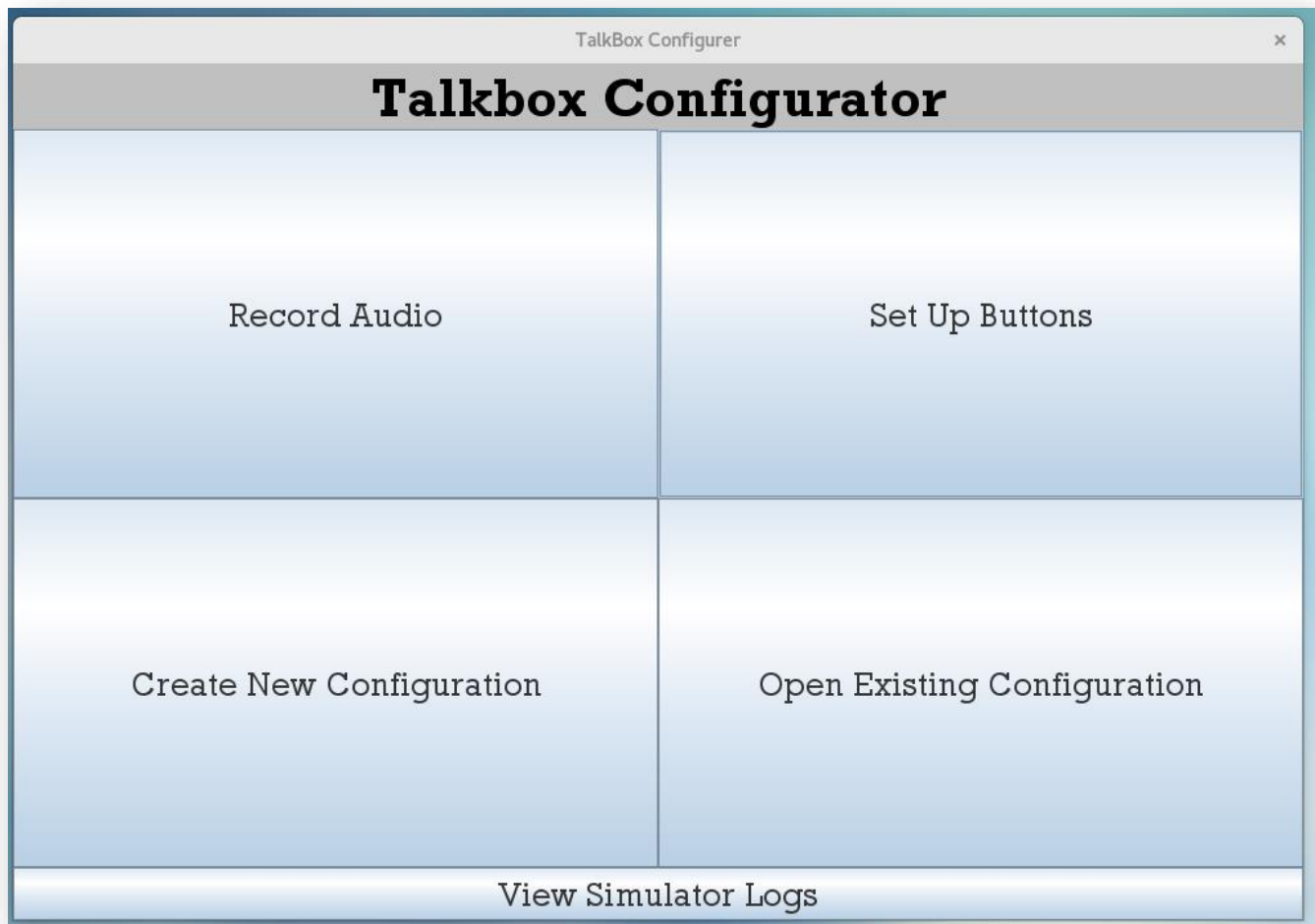

*All five buttons fully lit up.*

Talkbox Project  
EECS 2311  
Alberto M., Rohan T., Jordan M.


By default, there are six buttons, but the user can add as many buttons as they want through the “addition” of rows. The user can also delete rows that are not needed or not to their liking. In this panel the user can also record audio and add it straight to a button. Adding an image to the button is also a feature which is very user friendly. The user can select the orientation of the image (square, vertical or horizontal).

*Visual of the setup prompt for each button. Images and sounds can be added instantly.*




The 'Setup Button' dialog box is a window with a title bar containing a close button (X). It contains the following elements:

- Button Text:** A text input field with the value 'Button' and an **Add Emoji** button to its right.
- Sound Controls:** A section with a **Play Sound** button on the left and two buttons on the right: **Select Sound** and **Record Sound**, both in red text.
- Sound Path:** A text field displaying '(none)'.
- Select Color:** A large light blue button.
- Current Color:** A text field.
- Select Image:** A large light blue button.
- ImagePath:** A text field displaying '(none)'.
- Confirm Setup:** A large light blue button at the bottom.

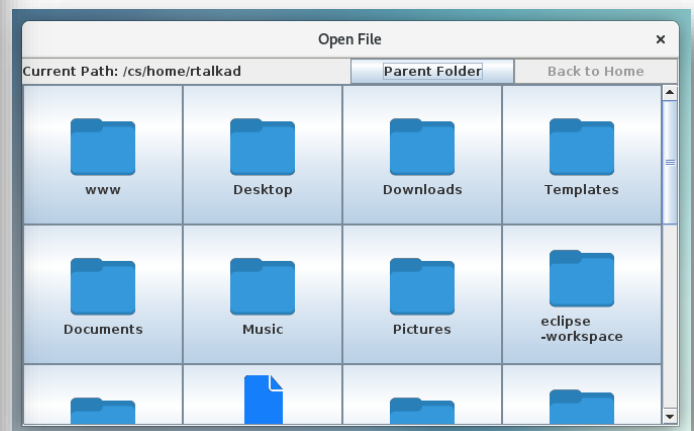

Talkbox Project  
EECS 2311  
Alberto M., Rohan T., Jordan M.


*This is the emoji panel. They can be dynamically added to buttons if one want to change its expression.*

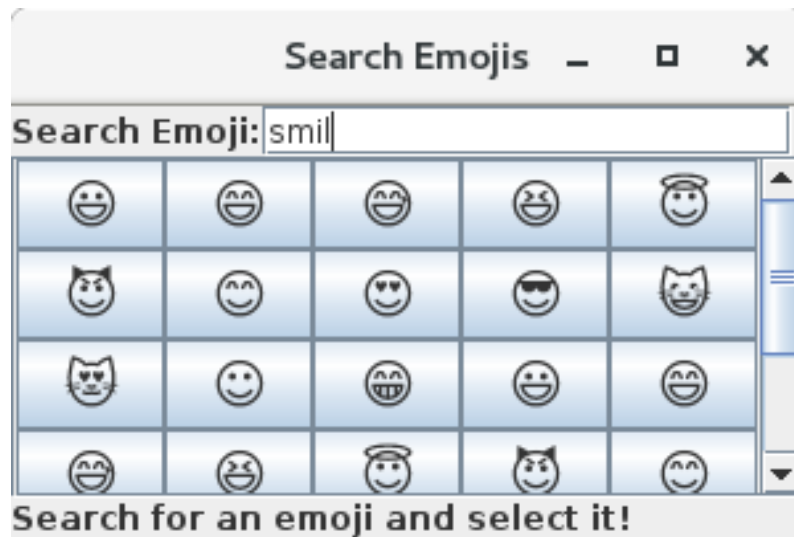

*Of course, the user can change the button color as they see fit.*



*The complete panel with a sample selection of images, sounds, and colors.*

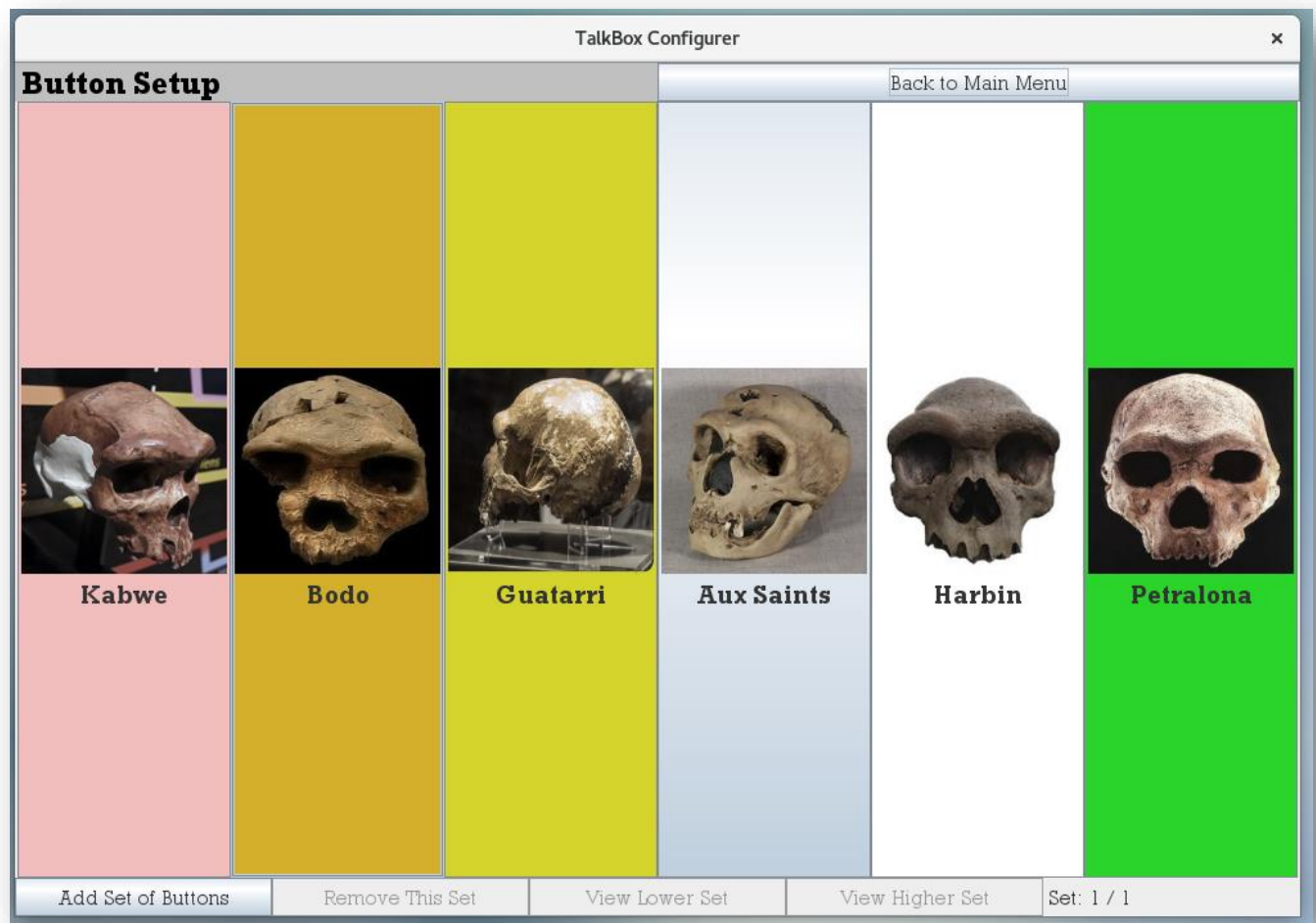

### 3 Viewing the Simulator Logs

Perhaps you, as a caregiver, are interested in the use of the simulator and whether certain buttons are being used. Maybe, there need to be different buttons for the same phrase. We record a every single user

Alberto M., Rohan T., Jordan M.

TalkBox Configurer

# Talkbox Configurator

Record Audio

ns

Create New Configuration

Open Existing Configuration

View Simulator Logs

**TBCLog**

**Menu**

- Apr 01, 2019 5:08:45 PM main.java.Talkbox.log.LogController logMess
- INFO: Set configuration for Guatarri
- Apr 01, 2019 5:08:45 PM main.java.Talkbox.log.LogController logMess
- INFO: Set configuration for Aux Saints
- Apr 01, 2019 5:08:45 PM main.java.Talkbox.log.LogController logMess
- INFO: Set configuration for Harbin
- Apr 01, 2019 5:08:45 PM main.java.Talkbox.log.LogController logMess
- INFO: Set configuration for Petralona
- Apr 01, 2019 5:08:45 PM main.java.Talkbox.log.LogController logMess
- INFO: Config added at /cs/home/rtalkad/Desktop/TalkboxData
- Apr 01, 2019 5:08:56 PM main.java.Talkbox.log.LogController logMess
- INFO: Select configuration button pressed.
- Apr 01, 2019 5:08:59 PM main.java.Talkbox.log.LogController logMess
- INFO: Simulator opened.
- Apr 01, 2019 5:09:00 PM main.java.Talkbox.log.LogController logMess



Talkbox Project

EECS 2311

Alberto M., Rohan T., Jordan M.